



IsoFORM® Back Braces

Comfortable compression for back stability and support

Designed for superior lower back protection and stabilization, the new DonJoy® IsoFORM® back braces help provide support and pain relief with customization that's as unique as your patients. They can aid with a comprehensive range of chronic and acute indications; supplement support with hot/cold packs for additional relief.

Built with Customizable Support, Targeted Compression & Patient Comfort in mind



· Adjustable Belt Wings – accommodates different shapes for a custom fit



· Paneling – provides durable rigid support



 Dual Pulley System – provides targeted compression with independent top and bottom compression



• Flexible Pull Tabs – for ease of use with adjustable lace length feature



• Donning Mitt – assists with ease of application



Semi-Universal – comfortably fits waist sizes from 28-60"



• Hexagonal Spacer Mesh – durable and breathable materials for comfort



• Optional Hot/Cold Pack – for additional relief



• EVA lined Panels – provides comfortable support



 Enclosed Lacing panel - protects lacing system from unwanted snags or tears



IsoFORM® - SIO

Flexible, low profile belt with single pulley system provides compression to reduce SI joint pain & dysfunction. Designed with silicone grip banding to help prevent brace from migrating. HCPC: L0621

Indications: suitable for people with sacroiliac pain, dysfunction, lower back pain, sprains/strains



IsoFORM - LO/LO+

The IsoFORM LO provides sagittal control and support from L1-L5. With optional Hot Cold Pack for additional relief. HCPC: L0626 (CST) / L0641 (OTS)

IsoFORM LO+ similar design as the LO, plus a Removable Semi-Rigid Anterior Panel for additional trunk support. HCPC: L0627 (CST) / L0642 (OTS)

Indications: Designed to support and stabilize the lower back and may be suitable for people with; acute and chronic low back pain, lower back sprains or strains, osteoporosis, spondylolisthesis



IsoFORM - LSO/LSO+

The IsoFORM LSO with Posterior Panel and removable Anterior Panel provides sagittal control and support from T9-S1. HCPC: L0631(CST) / L0648(OTS)

IsoFORM LSO+ Similar design as the LSO, plus two removable semi-rigid lateral extensions to help reduce lateral movement. HCPC: L0637(CST)/ L0650 (OTS) *Indications: Designed to support and stabilize the lower back and may be suitable for people with acute and chronic low back pain, lower back sprains or strains, osteoporosis, spondylolisthesis*

IsoFORM® - Postural Extension

The IsoFORM Postural Extension – provides comfortable thoracic region stability with a telescoping malleable aluminum posterior support & adjustable shoulder straps for a customized fit

- · Removable semi-rigid anterior panel for additional trunk support
- · Sternal Cross Strap secures the shoulder straps for a comfortable fit
- Adjustable backpack-style straps for easy donning and doffing HCPC: L0456 (CST) /L0457 (OTS)

Indications: suitable for people with thoracic spine conditions such as kyphosis, osteoporosis, acute pain, mechanical back pain, compression fractures, strain and mild scoliosis.











- A. Adjustable Straps
 - Backpack style for easy donning and doffing
- B. Removable Semi-Rigid Anterior Panel For additional trunk support
- C. Sternal Cross-Strap
 For a comfortable fit
- D. Adjustable Belt Wings

Accommodates for different shapes and allows clinician to customize the fit for a wide range of patients



- E. Telescoping Aluminum Posterior Support
 - Height adjustable and malleable to help tailor to the patient
- F. Enclosed Lacing Panel
 - Protects lacing system from unwanted snags or tears
- G. Dual Pulley System
 - With independent top and bottom compression
- H. Pull Tabs

For comfortable donning and doffing



IsoFORM® - TLSO

The IsoFORM TLSO – provides trunk restriction and stability from T3-S1 in three planes of motion with its overlapping panel system; Anterior, Posterior and Lateral Extension panels and a telescoping aluminum sternal bar that is adjustable and malleable for a custom fit.

- · Anterior panel helps restrict flexion and trunk rotation and supports T3-S1
- · Posterior panel rigid paneling helps restrict extension, trunk rotation and supports T3-S1
- · Lateral Extensions Removable semi-rigid supports help reduce lateral movement
- $\cdot\;$ Belt Wings with removable semi-rigid paneling encompasses the torso for customizable support
- · Easy to apply and remove with Pivoting Quick-release buckles
- · Over-the shoulder or Underarm Strap configuration options

HCPC: L0464

Indications: suitable for people with compression fractures of the thoracic spine, osteoporosis, Thoracic mechanical back pain and Kyphosis







A. Pivoting Quick Release Buckles

For easy donning and doffing

B. Aluminum Sternal Bar

Telescoping bar can be adjusted to desired height & malleable for custom fitting to desired angle

C. Semi-Rigid Belt Wings

Encompasses the torso for added support (removable plastic panel for customizable support)

D. Lateral Extensions

Semi-rigid supports to help reduce lateral movement (Removable for customizable support)

E. Anterior Panel

Helps restrict forward flexion and trunk rotation

F. 2 Strapping Configurations

Optional over the shoulder or underarm strap configuration

G. Dual Pulley System

Dual Pulley System with independent top and bottom compression

H. Enclosed Lacing Panel

Protects lacing system from unwanted snags or tears


I. Pull tabs

For comfortable donning and doffing

J. Posterior Panel

Plastic posterior panel helps restrict motion

Range of Spinal Support

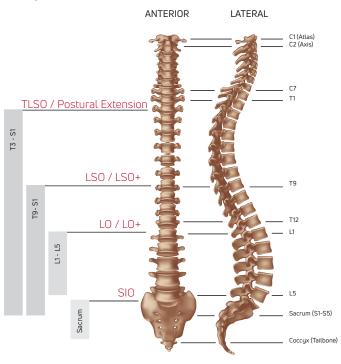

| DESCRIPTION | IsoFORM® SIO | IsoFORM LO / LO+ | IsoFORM LSO / LSO+ | IsoFORM Postural Extension | IsoFORM TLSO |
|-------------|-----------------|---------------------|-----------------------|----------------------------------|-----------------|
| S/M | 28″-48″ | 28"-48" | 28″-48″ | 28″-48″ | 28″-48″ |
| | (71cm-122cm) | (71cm-122cm) | (71cm-122cm) | (71cm-122cm) | (71cm-122cm) |
| L/XL | 48" – 60" | 48" – 60" | 48" – 60" | 48" – 60" | 48" – 60" |
| | (122cm – 152cm) | (122cm – 152cm) | (122cm – 152cm) | (122cm – 152cm) | (122cm – 152cm) |

Ordering Information

| cription | Size |
|------------------------|--|
| ORM® SIO | S/M |
| ORM SIO | L/XL |
| ORM LO | S/M |
| ORM LO | L/XL |
| ORM LO+ | S/M |
| ORM LO+ | L/XL |
| ORM LSO | S/M |
| ORM LSO | L/XL |
| ORM LSO+ | S/M |
| ORM LSO+ | L/XL |
| ORM Postural Extension | S/M |
| ORM Postural Extension | L/XL |
| ORM TLSO | S/M |
| ORM TLSO | L/XL |
| Cold Pack Accessory | Can be used on LO through TLSO product offerings |
| nal Support Accessory | XS/S |
| nal Support Accessory | XL/XXL |
| | ORM® SIO ORM SIO ORM LO ORM LO ORM LO+ ORM LO+ ORM LSO ORM LSO+ ORM LSO+ ORM LSO+ ORM Postural Extension ORM Postural Extension ORM TLSO Cold Pack Accessory mal Support Accessory |



