

enovis™

DYNACLIP QUATTRO™


ACTIVE, ADAPTIVE HEALING FOR MIDFOOT FUSION

The staple of your OR, DynaClip Quattro™ combines the strength of a plate with the active compression and ease of use of a staple.




STAPLE SYSTEM



The DynaClip Quattro Bone Fixation System excels in performance, fit, and reliability. Designed for the toughest fusion applications, the 4-leg configuration with an extra broad bridge provides additional compression and rotational stability, representing the **ultimate choice** in robust NiTiNOL staples for midfoot fusion applications.

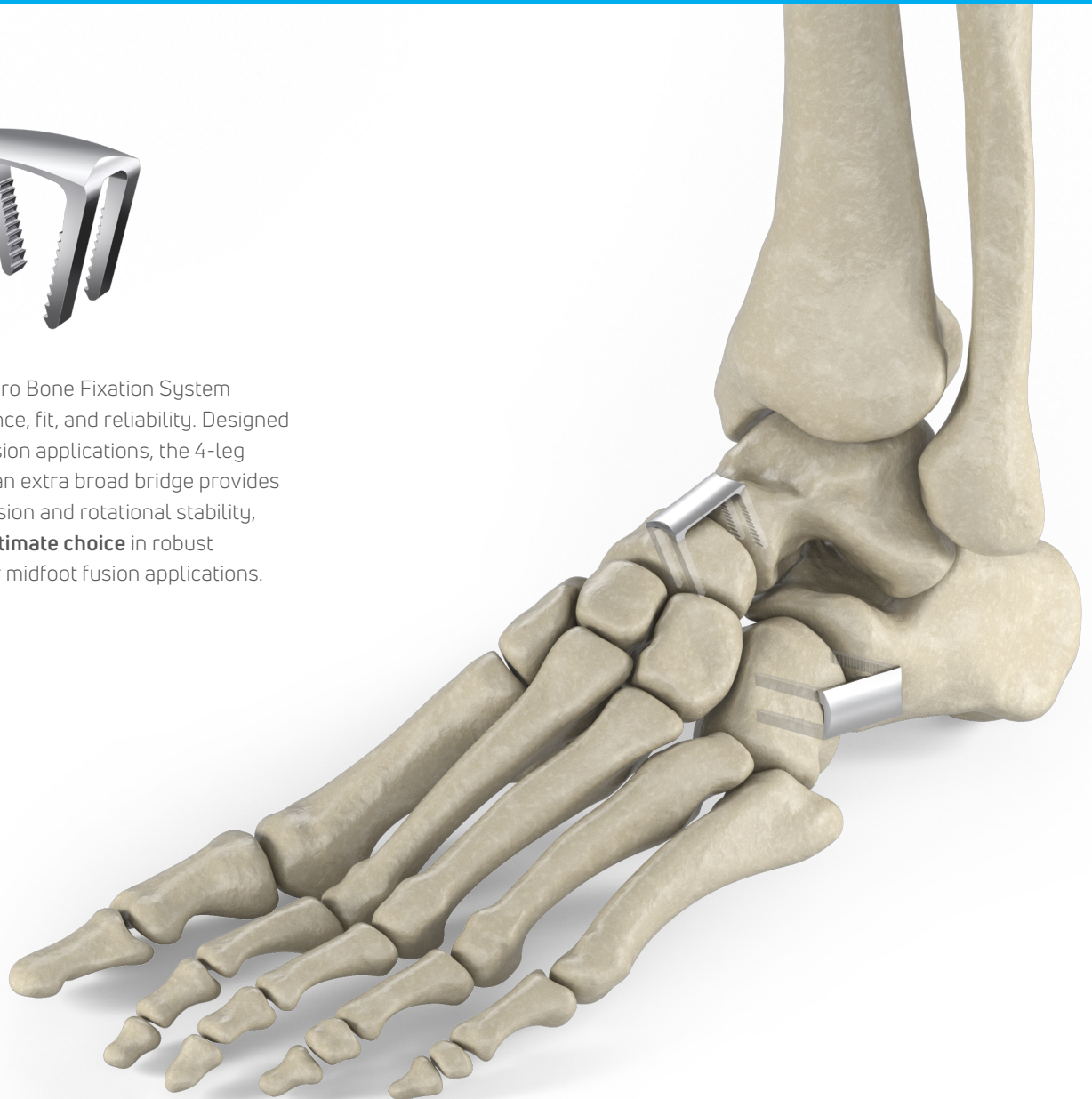

FEATURES & BENEFITS

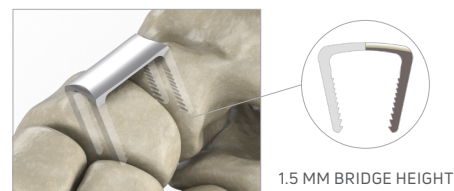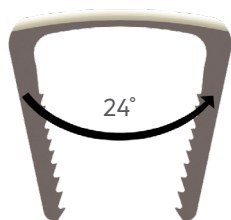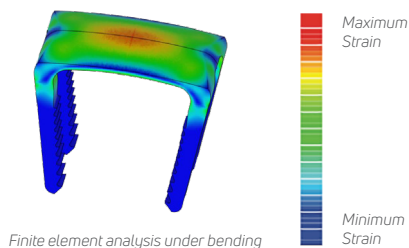

ENHANCED DURABILITY

DynaClip Quattro's design places stresses across the bridge, away from the legs, improving resistance to failure post-surgery.¹

MAXIMIZED COMPRESSION

Increased bridge curvature allows for compression generation through bridge, not just legs.¹

LOW PROFILE

Low-profile bridge eliminates any prominence or patient irritation.

EFFICIENT, ERGONOMIC INSTRUMENTATION

LOW-PROFILE DISPOSABLE INSERTER

DynaClip Quattro comes pre-loaded on a disposable Inserter that empowers:

- Simple, intuitive deployment.
- Accurate leg positioning every time.
- Fewer instruments. Just insert, deploy, and tamp using the disposable inserter.



A COMPREHENSIVE SYSTEM

The DynaClip family offers many options for a wide range of indications based on surgeon preference and patient anatomy.



References:
1. Data on File, MedShape, Inc., 2022

enovis

T 800.495.2919 F 877.778.3864

Medshape, Inc.
1575 Northside Drive NW | Suite 440 | Atlanta, GA 30318 | U.S.A.
enovis.com/foot-and-ankle

Copyright © 2022 DJO Foot and Ankle
MK-10273 Rev 00

DYNACLIP QUATTRO IMPLANTS & SINGLE-USE INSTRUMENTS

PART #	DESCRIPTION
3000-11-181616	DYNACLIP QUATTRO, 18 MM X 16 MM X 16 MM
3000-11-221616	DYNACLIP QUATTRO, 22 MM X 16 MM X 16 MM
3000-11-222020	DYNACLIP QUATTRO, 22 MM X 20 MM X 20 MM

Individual results may vary. Neither Medshape, Inc. nor any of the Enovis companies dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.