enovis

Forte (n.): Strength, excels in a specific area

DYNACLIP FORTE®

BONE FIXATION SYSTEM




The DynaClip Forte® Bone Fixation System

excels in performance, comfort, and reliability. Designed for the most demanding weight-bearing applications, the 4-leg configuration with a broad bridge provides additional compression and rotational stability, representing the ultimate in NiTiNOL staples.



FEATURES & BENEFITS

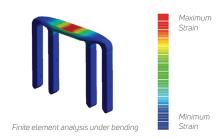

ENHANCED DURABILITY

DynaClip Forte's patented design places stresses across the bridge and away from the legs, improving resistence to failure post-surgery.¹



MAXIMIZED COMPRESSION

Increased bridge curvature allows for compression generation through bridge, not just legs.¹





BRIDGE HEIGHT

LOW PROFILE

Low-profile, tapered bridge eliminates any prominence or patient irritation.

EFFICIENT INSTRUMENTATION

DISPOSABLE PROCEDURE PACK

All instrumentation comes packaged in one sterile kit, designed to eliminate waste and improve efficiency in the OR. No need to open multiple packs.

- **1.** 3mm Drill (1)
- 2. 3mm Locator Pins (4)
- 3. Universal Drill Guide





LOW-PROFILE DISPOSABLE INSERTER

DynaClip Forte® comes pre-loaded on a disposable Inserter that empowers:

- · Simple, intuitive deployment.
- · Accurate leg positioning every time.
- No adjusting needed in the OR. Simply insert and deploy.



INSERT DYNACLIP FORTE® PULL UP

PULL UP & SLIDE TO RELEASE

TAMP USING INSERTER

A COMPREHENSIVE SYSTEM

The DynaClip family offers many options for a wide range of indications based on surgeon preference and patient anatomy.







DYNACLIP FORTE IMPLANTS & SINGLE-USE INSTRUMENTS

 DESCRIPTION
 PART #

 DYNACLIP, 4-LEG INLINE, 26 MM x 16 MM
 3000-06-261616

 DYNACLIP, 4-LEG INLINE, 26 MM x 20 MM
 3000-06-262020

 DYNACLIP, 4-LEG INLINE, 30 MM x 16 MM
 3000-06-301616

 DYNACLIP, 4-LEG INLINE, 30 MM x 20 MM
 3000-06-302020

 DISPOSABLE PROCEDURE PACK, STERILE
 3000-04-000

References: 1. Data on File, MedShape, Inc., 2023

enovis

T 800.495.2919 F 877.778.3864

Medshape, Inc.
1575 Northside Drive NW | Suite 440 | Atlanta, GA 30318 | U.S.A. enovis.com/foot-and-ankle

Copyright © 2023 MedShape, Inc. MK-10189 Rev 01 Individual results may vary. Neither Medshape, Inc. nor any of the Enovis companies dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.